Evolve Scope Drying Cabinet

Storage Solutions for Healthcare






Designed to support a hospital's efforts to decrease infection and improve patient safety

Effective endoscope drying is essential to the prevention of bacterial transmission and hospital-acquired infections. Endoscope storage guidelines from leading organizations such as SGNA, AORN, and AAMI cite the importance of keeping endoscope channels dry between reprocessing and use. Any presence of moisture can lead to microbial contamination. Our scope cabinets provide ventilated, clean, cushioned, and dust-free environments for storing endoscopes, but our Evolve scope drying cabinet includes an automated system that provides consistent and continuous delivery of HEPA-filtered, pressurized air into all endoscope channels to keep scopes dry and avoid reprocessing.

The easy-to-operate drying system features an automated, 10-minute drying cycle per scope and continuous venting cycle for overnight storage. It has a self-adjusting, high-pressure regulation and automatic shut-off and offers a HEPA-filtered 99.97% retention rate of particles less than 0.3µm. The hands-free unit works with major endoscope brands.

The cabinet is constructed of proprietary AireCore, a durable, nonporous, and bacteria- and moisture-resistant material that further supports a hospital's infection prevention goals. AireCore also offers customization capabilities and design flexibility. AireCore surface finishes include Glacier White, Brushed Aluminum, or specified laminate. Cabinets are also available in stainless steel.

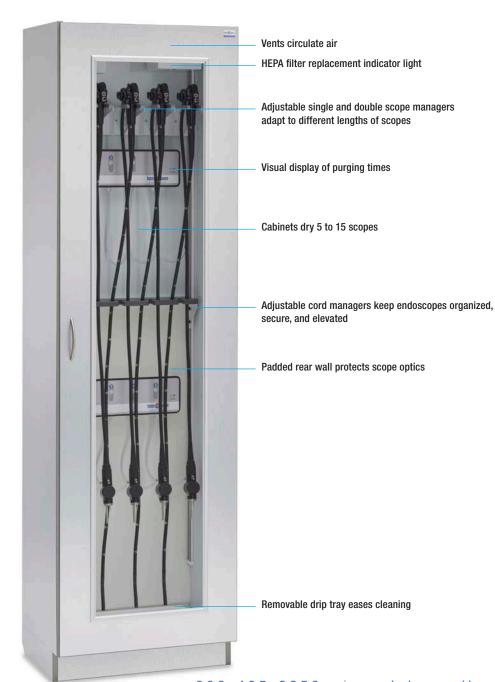







Prepare your endoscopy department now for future changes to compliance standards with a scalable solution that addresses your department's standards and budgets. Our scope drying unit can be ordered separately and field-installed in a new or previously purchased Evolve scope cabinet with AireCore or stainless steel construction.

The drying system requires power input of 90-250 VAC and 50-60 Hz; maximum power output is 12 VDC, 7.5 amps, and 90 watts. The cabinet includes a 6-foot power cord. Disposable connectors must be purchased separately.



Exterior (interior) dimensions:

26.75"w x 19.25"d x 92"h (25.25"w x 17"d x 86.5"h)
92"h cabinet with 6"h toe kick is 94"h

Scope capacity:

Single scope manager holds 5 scopes Double scope manager holds 9 scopes 83" maximum scope length

Interior Configuration:

AireCore:

SEA1927SGLSF, SEA1927SGRSF; with double scope holder: SEA1927SGLDF, SEA1927SGRDF

AireCore with Brushed Aluminum:

SEB1927SGLSF, SEB1927SGRSF; with double scope holder: SEB1927SGLDF, SEB1927SGRDF

Stainless Steel:

SES1927SGLSF, SES1927SGRSF;

with double scope holder: SES1927SGLDF, SES1927SGRDF







Exterior (interior) dimensions:

36"w x 19.25"d x 92"h (34.5"w x 17"d x 86.5"h) 92"h cabinet with 6"h toe kick is 94"h

Scope capacity:

Single scope manager holds 8 scopes

Double scope manager holds 15 scopes

83" maximum scope length

Interior Configuration:

AireCore:

SEA1936SGSF; with double scope holder: SEA1936SGDF

AireCore with Brushed Aluminum:

SEB1936SGSF; with double scope holder: SEB1936SGDF

Stainless Steel:

SES1936SGSF; with double scope holder: SES1936SGDF

The cabinet is warranted for 5 years; the drying unit is warranted for 1 year.

UL recognized, CE marked, meets Medical Safety Standard 60601-1, 4th edition.





Solaire Medical specializes in storage solutions for all areas of your healthcare facilities. We offer the InnerSpace® portfolio of products, a cohesive system of carts, cabinets, and modular casework that improves staff and facility efficiency and an organization's bottom line. Our InnerSpace portfolio is supplemented with the Quick line of open storage products.

When you purchase InnerSpace by Solaire Medical products, you're buying products made in the U.S.A. and from a direct manufacturer that delivers on-time and error-free. The portfolio's five-year

warranty is our promise to you that our products and service are of the highest quality.

The Quick family of open storage products have a one-year warranty.

Contact our service team with your requests or to connect with a sales representative. We're here to help you specify the best storage products for your needs, whether a standard option or customized solution. We offer complimentary storage analysis, design services, CAD support, and logistics and installation assistance.

