

Roam Supply Carts

Storage Solutions for Healthcare



Smart design. Flexible design. Durable design.

Roam's efficient exterior footprint and next-to-no wasted interior space result in more storage capacity per square inch. Roam will perform day in and day out and look good doing it.







Roam 1 counter-height cart

General storage cart configuration includes 1 single-compartment shelf, 2 shelf dividers, 1 4" basket, 1 8" basket, and label holders (SR1CGGS)



Roam 2 counter-height cart

General storage cart configuration includes 1 single-compartment shelf, 2 shelf dividers, 1 4" basket, 1 8" basket, 4 4" trays with dividers, 1 8" tray with dividers, and label holders (SR2CGGS)

Smart design

The smart design of the FlexCell system and interior components maximize Roam's storage capacity, whether you are outfitting a counter-height or full-height cart with one, two, three, or four compartments.

Glass doors are designed with a position hinge that holds the door open to exactly where you want it to be. A magnetic catch securely closes and keeps the door from opening in transit.

The roll-top door track is located on the outside of the cart, so interior storage space is not compromised.



Roam 1 cart

General storage cart configuration includes 1 single-compartment shelf, 2 shelf dividers, 2 2" trays with dividers, 3 4" trays with dividers, 2 4" baskets, 1 8" basket, 1 8" tray with dividers, and label holders (SR1RGS)



Roam 2 cart

ENT/Neuro cart configuration includes 4 single column shelves, 8 shelf dividers, 4 2" trays with dividers, 6 4" trays with dividers, 2 4" baskets, 2 8" baskets, and label holders (SR2GNU)

Durable design

Roam takes the wear and tear that mobile supply carts get. A metal base and top add stability to the case. Plastic ABS interior components won't wear down FlexCell tracks, even with frequent pulling out and pushing in of baskets and trays.


Two roll-top doors on our Roam 3 and Roam 4 make it easier to open and close and also eliminate damage common with wider roll-top doors. Both roll-top and hinged doors operate smoothly.

Even our widest Roam carts are still mobile carts, and they need to be manageable and durable. A lightweight aluminum composite frame aids maneuverability, and soft bumpers protect carts, walls, and doors against damage.



Roam 3 cart

Boxed catheter and general storage cart configuration includes 2 single-compartment shelves, 2 4" trays with dividers, 1 8" tray with dividers, 2 8" baskets, 1 double-compartment shelf, 9 shelf dividers, 1 boxed catheter module, and label holders (SR3RBG)

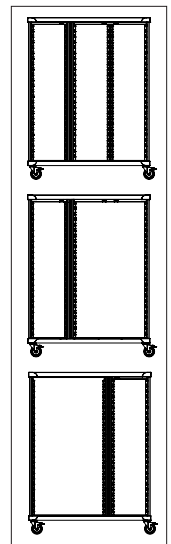

Compartment configurations

Flexible design

The FlexCell design and interior components create a modular system that can be configured and reconfigured as your needs and supplies change. Center columns can be retrofit to carts.

Tracks on the FlexCell system provide two tilt options for trays and baskets: a moderate tilt for most applications or a gravity-feed tilt for ER or OR applications where quick access is paramount.



Roam 4 cart


General storage cart configuration includes 8 single-compartment shelves, 16 shelf dividers, 4 2" trays with dividers, 16 4" trays with dividers, 4 4" baskets, and 4 8" baskets (SR4GGS)

Clean design

Supply carts must be highly functional, but that doesn't mean they can't look good. Concealed hinges and a flush door design add up to a better aesthetic.

Nonporous surfaces, tight seams, and few exposed fasteners support infection prevention goals and reduce dust and dirt traps.

The FlexCell system and interior components are designed for ease of cleaning. Roam carts can be cleaned with agents commonly used in hospitals. Aluminum composite surfaces won't rust or scratch.

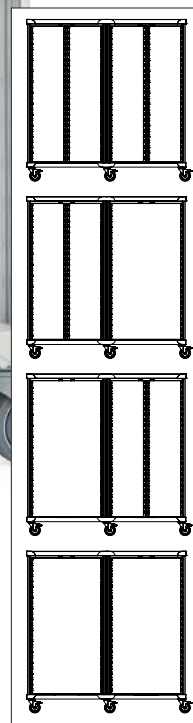

Compartment configurations



Product shown: Roam 2 cart with center column and accessories (SR2G) in Wilsonart 5th Avenue Elm 7966



Position hinge keeps door open to the position user places it



Electronic lock and keyed lock options;
no lock is standard



Easy-grip stainless steel push/pull handle



Exterior side-mounted accessories add flexibility
to Roam carts



5" lockable medical-grade casters



Product shown: Roam 2 cart with accessories (SR2R)

Maximum Versatility

Maximize interior storage with FlexCell and a modular system of trays, baskets, and dividers.



2"-, 4"-, and 8"-high trays let you match the size of supplies to the size of the tray and reduce unused vertical space.



Stackable trays make transporting restocked trays efficient.



Powder-coated standard metal shelf spans one compartment and holds up to three cath managers®; multipurpose shelf comes in a single-compartment shelf that holds up to three cath managers and a double-compartment shelf that holds up to seven cath managers.



Labels clip on trays and baskets help organize restocking and reduce time wasted searching for supplies.



Dividers are adjustable and removable so you can organize and reorganize tray contents. Clips keep dividers securely in place.



Safety stops on the FlexCell tracks and backs of trays and baskets prevent them from sliding out of cart under normal use.



4"- and 8"-high clear baskets help you identify contents quickly. Perforated bottoms prevent dust build up.



Trays and baskets have a moderate and gravity feed tilt and are designed with straight, not sloped sides, resulting in a higher percentage of usable space and more capacity. Our 4"-high tray, compared with a sloped 4"-high tray, provides up to 38% more storage volume.



Product shown: Roam 2 long and short catheter cart with optional embolectomy rack (SR2GLSC))

Maximum Functionality

Customize your Roam cart with an offering of accessories mounted to the cart's exterior.



9-hook cath managers mount on multipurpose shelves. Each hook holds 10 to 15 catheter packages. Includes product labeling. (S9CM)



Side-mounted rail is available on standard and counter-height carts.



Pull-out suture modules include shelves and dividers designed to hold and access suture boxes and other small supplies. Each shelf holds 8 2.25" deep suture boxes for a total of 48. (SPSM)



Full-height carts are available with sloped tops. All carts are available in Brushed Aluminum laminate.



Easy-to-clean aluminum slat wall holds a variety of accessories and is available on all Roam carts.

Roam Product Line

Two heights, four widths, and a straightforward kit of components and procedure-specific accessories let you create exactly the cart you need.



Roam 1 counter-height cart

The Roam 1 counter-height cart has a glass door that can be hinged on the right or left side. The cart features a metal work surface.

Exterior dimensions:

20"w x 28.75"d x 44.25"h

Weight: 95 lbs

Product Number: SR1CGL, SR1CGR

Roam 2 counter-height cart

The Roam 2 counter-height cart has two lockable tempered glass doors. The cart can be ordered with or without a center column.

Exterior dimensions:

37.5"w x 28.75"d x 44.25"h

Weight: 145 lbs*

Product Number: SR2CG

Roam 1 cart

The Roam 1 cart comes with one lockable counterbalanced roll-top door or one lockable tempered glass door that can be hinged on the right or left side. A 5" sloped top is also available on carts with glass doors.

Exterior dimensions:

20"w x 28.75"d x 75.25"h (glass door)

20"w x 28.75"d x 81"h (roll-top door)

Weight: 120 lbs

Product Number:

SR1GR, SR1GL (glass door)

SR1R (roll-top door)

Roam 2 cart

The Roam 2 cart has one lockable counterbalanced roll-top door or two lockable tempered glass doors. The cart can be ordered with or without a center column.

Exterior dimensions:

37.5"w x 28.75"d x 75.25"h (glass doors)

37.5"w x 28.75"d x 81"h (roll-top door)

Weight: 193 lbs*

Product Number:

SR2G (glass doors)

SR2R (roll-top door)

Roam 3 cart

The Roam 3 cart has two lockable counterbalanced roll-top doors or a single and double set of lockable tempered glass doors. The cart can be ordered with a single compartment on the left and two center columns, with a right center column, or with a left center column.

Exterior dimensions:

57"w x 28.75"d x 75.25"h (glass doors)

57"w x 28.75"d x 81"h (roll-top doors)

Weight:

293 lbs (glass doors)*

288 lbs (roll-top doors)*

Product Number:

SR3G (glass doors)

SR3R (roll-top doors)

Roam 4 cart

The Roam 4 cart has two lockable counterbalanced roll-top doors or two sets of lockable tempered glass doors. The cart can be ordered with three center columns, two center columns on the right side, two center columns on the left side, or one center column.

Exterior dimensions:

74"w x 28.75"d x 75.25"h (glass doors)

74"w x 28.75"d x 81"h (roll-top doors)

Weight:

376 lbs (glass doors)*

371 lbs (roll-top doors)*

Product Number:

SR4G (glass doors)

SR4R (roll-top doors)


Soft White is the standard paint finish on Roam carts. Five additional paint colors and a wide range of laminates are available.



*Weights shown are for carts with center columns. Counter-height center column is 10 lbs; full-height center column is 23 lbs.




2" Tray with Dividers (STD2)

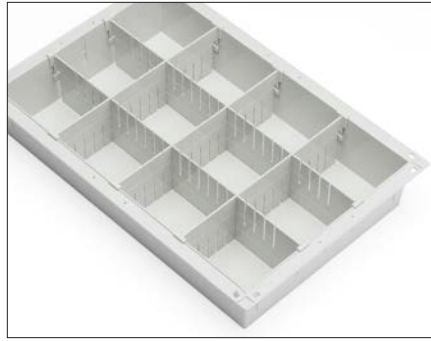

4" Tray with Dividers (STD4)



8" Tray with Dividers (STD8)

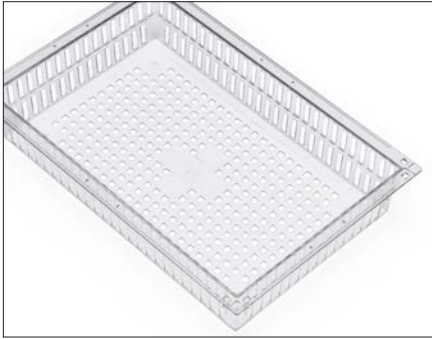


4" Basket (SB4)



8" Basket (SB8)



4" and 8" Basket Dividers (SBD4, SBD8)

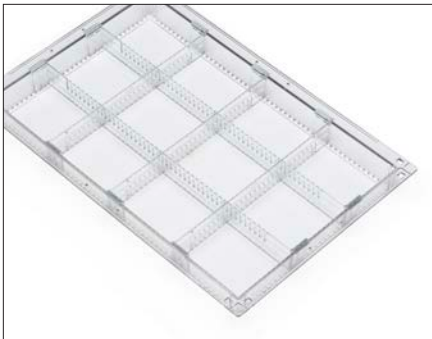

2" Clear Tray with Dividers (SCT2)



Label Holders (SLH35)



Metal Shelf (S27S1)



Multipurpose Shelf (S27M1) with 3-Cath Manager (S9CM)



Lockable Drawer (single compartment, S27D1)



Multipurpose Shelf (S27M1) with 2 Dividers (S27DI)



Multipurpose Shelf (S27M2) with 7-Cath Manager (S9CM)



Lockable Drawer (double compartment, S27D2)



Multipurpose Shelf (S27M2) with 5 Dividers (S27DI)

Preconfigured Carts

Here is a sampling of Roam carts preconfigured to address the specific storage needs of departments throughout a healthcare facility.

Cath Labs and Radiology Departments



Roam 2 Catheter Cart

Includes 1 double-compartment shelf, 5 shelf dividers, and 7 9-hook cath managers. Glass (SR2GCC) or roll-top (SR2RCC) door



Roam 2 Boxed Catheter Cart

Includes 1 double-compartment shelf, 5 shelf dividers, and a pull-out module with adjustable dividers to hold catheter boxes. Pull-out module is divided to prevent tipping. Glass (SR2GBC) or roll-top door (SR2RBC)



Roam 3 Catheter Cart

Includes 1 single-compartment shelf, 1 double-compartment shelf, 7 shelf dividers, and 10 9-hook cath managers. Glass (SR3GCC) or roll-top (SR3RCC) doors



Roam 3 3-Manager Catheter and Supply Cart

Includes 4 single-compartment shelves, 8 shelf dividers, 3 9-hook cath managers, 6 4" trays with dividers, 4 8" trays with dividers, 2 8" baskets, and label holders. Glass (SR3G3CG) or roll-top (SR3R3CG) doors

Surgery Departments



Roam 2 Laproscopy Cart

Includes 4 single-compartment shelves, 8 shelf dividers, 8 4" trays with dividers, 2 8" baskets, and label holders. Shown in brushed aluminum. Glass (SR2GLAP) or roll-top (SR2RLAP) door



Roam 2 Open Heart Cart

Includes 3 single-compartment shelves, 6 shelf dividers, 1 2" tray with dividers, 4 4" trays with dividers, 1 8" tray with dividers, 2 4" baskets, 1 8" basket, 1 pull-out suture module, and label holders. Shown in brushed aluminum. Glass (SR2GOH) or roll-top (SR2ROH) door



Roam 3 Vascular Cart

Includes 6 single-compartment shelves, 12 shelf dividers, 2 2" trays with dividers, 8 4" trays with dividers, 4 8" trays with dividers, 1 pull-out suture module, and label holders. Shown in brushed aluminum. Glass (SR3GVS) or roll-top (SR3RVS) doors



Roam 3 Cysto Cart

Includes 7 single-compartment shelves, 14 shelf dividers, 3 9-hook cath managers, 1 2" tray with dividers, 8 4" trays with dividers, 3 8" trays with dividers, and label holders. Shown in brushed aluminum. Glass (SR3GCY) or roll-top (SR3RCY) doors

General Storage Areas



Roam 1 General Storage

Includes 1 single-compartment shelf, 2 shelf dividers, 2 2" trays with dividers, 3 4" trays with dividers, 2 4" baskets, 1 8" basket, 1 8" tray with dividers, and label holders. Glass (SR1GRGS) or roll-top (SR1RGS) door



Roam 2 Cart with 5 Shelves

Includes 5 shelves with 5 dividers each (25 total). Glass (SR2GDS) or roll-top (SR2RDS) door



Roam 2 Endovascular Cart

Includes 3 single-compartment shelves, 4 shelf dividers, 3 9-hook cath managers, 3 4" trays with dividers, 3 8" baskets, and label holders. Glass (SR2GEND) or roll-top (SR2REND) door



Roam 2 Cysto Cart

Includes 5 single-compartment shelves, 3 9-hook cath managers, 8 shelf dividers, 2 4" baskets, 1 8" basket, 1 8" tray with dividers, and label holders. Glass (SR2GCY) or roll-top (SR2RCY) door



Roam 3 7-Manager Catheter and General Storage Cart

Includes 1 single-compartment shelf, 1 double-compartment shelf, 2 shelf dividers, 7 9-hook cath managers, 5 4" trays with dividers, 1 8" tray with dividers, 1 4" basket, 1 8" basket, and label holders. Glass (SR3G7CG) or roll-top (SR3R7CG) doors



Roam 3 Pack/Drape Cart

Includes 4 single-compartment shelves and 4 double-compartment shelves. Glass (SR3GPD) or roll-top (SR3RPD) doors



Roam 3 Ortho/Cast Cart

Includes 6 single-compartment shelves, 12 shelf dividers, 4 4" baskets, 2 2" trays with dividers, 6 4" trays with dividers, 2 8" trays with dividers, 2 8" baskets, and label holders. Glass (SR3GOC) or roll-top (SR3ROC) doors




Roam 4 General Storage Cart

Includes 8 single-compartment shelves, 16 shelf dividers, 4 2" trays with dividers, 16 4" trays with dividers, 4 4" baskets, and 4 8" baskets. Glass (SR4GGS) or roll-top (SR4RGS) doors

Specialty Carts


In addition to our standard Roam offering, we offer two specialty carts that address unique procedure-specific equipment and supplies.

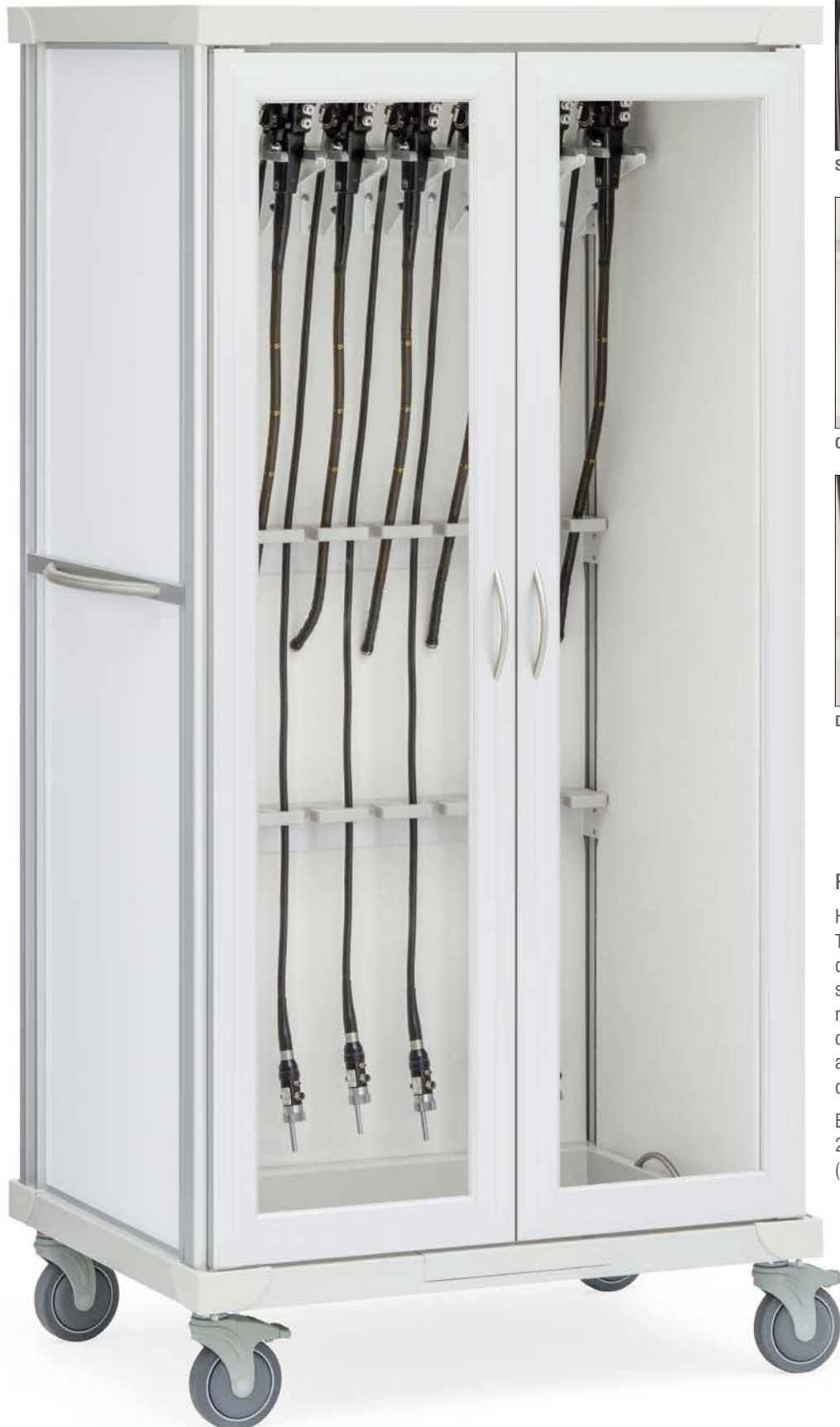

Scope Manager®



Cord Retainer

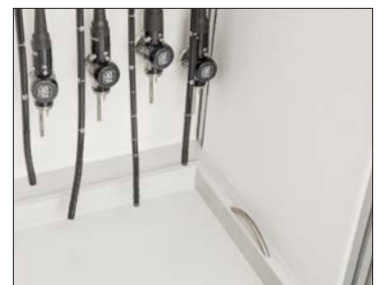

Drip Tray

Roam 2 Scope Cart

Holds 8 endoscopes up to 60" in length. The cart has lockable tempered glass doors, cushioned back panel to protect scope optics, ventilated bottom, and removable drip tray. Scope managers and cord retainers keep contents organized and elevated off the cart floor. Scope cart does not accept colonoscopes.

Exterior Dimensions:
28.75"d x 37.5"w x 75.25", 200 lbs
(SR2GS)



Roam 2 Long Catheter Cart

Cart is 80.5" high to accommodate catheters up to 64" in length. Comes with 7 9-hook cath managers permanently mounted to the top of the cart. Cath managers can be adjusted horizontally.

Exterior Dimensions: 28.75"d x 37.5"w x 80.5", 215 lbs (SR2GC)



Products shown: Front cover, left to right: Roam 2 EP cart (SR2GEP), Roam 2 counter-height general storage cart (SR2CGGS), Roam 3 arthroscopy cart (SR3RAR); Back cover, left to right: Roam 2 cart with accessories (SR2G) in Wilsonart Sunshine D499, Roam 1 PPE cart (SR1RPPE), Roam 2 counter-height cart (SR2CGNC)

InnerSpace

PORTFOLIO OF PRODUCTS

Solaire Medical specializes in storage solutions for all areas of your healthcare facilities. We offer the InnerSpace® portfolio of products, a cohesive system of carts, cabinets, and modular casework that improves staff and facility efficiency and an organization's bottom line.

When you purchase our products, you're buying products made in the U.S.A. and from a direct manufacturer that delivers on-time and error-free.

Contact our service team with your requests or to connect with a sales representative. We're here to help you specify the best storage products for your needs, whether a standard option or customized solution. We offer complimentary storage analysis, design services, CAD support, and logistics and installation assistance.

Call to schedule a free storage consultation.

solaire
MEDICAL

Solaire Medical | P.O. Box 2808 | Grand Rapids, Michigan 49501
888-435-2256 phone | 616-743-5944 fax | solairemedical.com



Roam
Supply Carts



Rover
Workstations



Tempo
Procedure Carts



Evolve
Cabinets



Evolve
Casework